

# Jenx Monkey Stander

Specifically designed for younger children, the Jenx Monkey is a gentle and fun introduction to standing, providing prone to upright positioning.

## Simple and highly adjustable standing system




**De-rotational thoracic & pelvis straps** with lateral pads that fasten in the centre and adjust independently from either side, enabling a mid-line position to be achieved and maintained

**Integral angle-adjustable tray** with both bowl and smooth tray infill for a wide variety of different activities

Polyurethane support surfaces **with antimicrobial agent**

**Quick, simple & easy** to use. Comes as a complete unit and very few adjustments are required



**Gentle, bright & friendly**

First diagnosis of a child's special needs is difficult for parents who have to learn a whole new vocabulary and deal with equipment and treatments that are new and bewildering. The Monkey has been carefully designed to offer the best in early intervention support in a friendly way.



Flexible growth: chest and hip boards can be independently altered and fixed to accommodate child.

Fully upholstered chest, hip and leg boards help make standing more comfortable, to increase standing tolerance.

Angle can be adjusted with the child in the frame and set to the correct angle with infinite adjustment.

## ACCESSORIES


### Tray Grab Bar

When 1 hand needs to fix to increase functionality  
5050-7201-000




### Activity Frame

Encourages head control & hand/eye coordination  
3615-7475-000

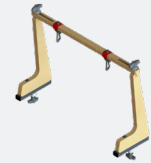

### Sandals

Size 1: 5356-6081-000  
Size 2 larger users or with splints: 5356-6082-000

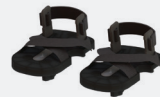

### Sandal Raising Blocks

Size 1: 5356-6091-000  
Size 2: 5356-6092-000



## CASE STUDY: CHILD A

**Condition:** Hydrocephalus, prematurity resulting in prolonged hospital stay from birth.

**Background:** Child A was born prematurely with complex needs resulting in a prolonged 18 month stay in hospital. Once home Child A was on home O<sup>2</sup>, PEG fed and presenting developmental delay. He has moderate head control, some trunk activity, but no independent sitting balance.

General and Emotional benefits (from trialling the Jenx Monkey): Child A is able to practice weight bearing and ensure better hip development and bone density. Plus the change of position gives Child A greater access to the world and interaction on that level. As Child A is very small even though nearly 2 years old, sandal raising blocks were needed to ensure the appropriate level of support and a good stand are achieved. The stander has also meant reduced physical support required to be provided by Mum, allowing her to play with Child A on better level, making both mum and Child A much happier.

Code	Age	Max User Weight	Angle Range	Base Size	Chest Height (Foot plate-Chest pad top)	Chest Pad Width	Hip Pad Width	Tray Size
5050-0001-000	9 mths-4 yrs	25 kg	20°-90°	88 x 57.5 cm	51-75 cm	66 x 50 x 60 cm	15-23.5 cm	36 d x 53 w cm