

Safe Surround Genie Bed

The Safe Surround Genie bed is for children and adults that can be found in dangerous positions due to movement during the night and need a safe environment for greater containment.

Superior Nursing Care for Children and Adults

Perspex panels provide free eye contact with carers, and the individual can take in what is going on around them

High quality craftsmanship: durable Beech timber with attractive aesthetics



For **carers in wheelchairs** Genie provides ample clearance of up to 74 cm

Sufficient clearance for **hoist transfers**

Standard Features



Adjustable lying surface including 4 section profile sleeping platform, Standard foam mattress, 4 locking castors.



Raise mattress level high to assist carers when tending to the client.



Lower mattress level to enable independent transfers.

SAFETY FEATURES FOR GENIE & PLUS BEDS


- Safe interior to prevent entrapment of head or limbs.
- Special door locking system prevents the doors being opened from the inside.
- Unique door hinge system that doesn't entrap fingers.
- High quality fittings and locking mechanisms.
- Power cut-off plug. When the hand switch is not activated (power disconnected state) all electrical components of the bed are disconnected from the mains and therefore no electricity is used. Electric and magnetic fields in the bed are prevented.
- Maintenance-free DC motors and manual safety switch.
- Emergency lowering in the case of power failure
- Special wooden surface coating that seals the surfaces for maximum hygiene.
- Sweat and saliva resistant wood treatment in accordance with DIN 71-3 Children's Toys Safety Act.

OPTIONS

Mattress Incontinence Cover
7956-0026-000



Pressure Relief Foam Mattress
7956-0011-051

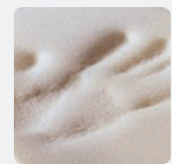

Internal Padding as Required. Velcro attachment.



Code	Door Height	Clearance Height (For Wheelchair)	Weight Limit	Inside Dimensions	Overall Height	Lying Height	Outside Dimensions
7959-0016-000	98 cm	74 cm	150 kg	200 x 90 cm	Min. 137 cm Max. 197 cm	Min. 42cm Max. 102 cm	223 x 101 cm

